

TIMERTAG



2019
PRODUCT RANGE

WWW.TIMERTAG.COM

TIMERTAG

GETTING CRITICAL TASKS DONE ON TIME CAN BE CHALLENGING MAKE IT EASY WITH TIMERTAG

Timertag's proven solutions use flashing lights to alert users about expiry and engage patients in their care. There's no sound, so there's no contribution to alarm fatigue. Available in single use or reusable formats, there's a Timertag solution to meet almost any need or budget.



LABELS

Use in a chart, at the bedside or on equipment to bring a task to attention.



WRISTBANDS

Perfect for alerting review of IV access placed outside a healthcare facility.



MINI TIMERTAG

A robust alternative to labels - ideal for equipment service alerts.



REUSABLE TIMERTAG

Our economical solution with unlimited use from 6 months to 5 years..

TABLE OF CONTENTS



NURSING

Emergency IV Access
IV Access Tubing Changes
Dressing Changes
IV Access Reviews
Oral Care
Urinary Catheter Reviews
Patient Repositioning
Create your own solution

PAGE

5
6
7
8 - 9
10
11
12
13



BIOMED

Mobile Equipment Service Alerts

15



PHARMA

Compounded Pharmaceutical Expiry Alerts

17



CSSD

Sterile Packs
Endoscope Storage
Medical Waste Containers

19
20
21



TIMERTAG

Independent Research
Awards
Environmental Statement
The Company
Our Team
Product Technical Specifications

22
23
24
25
26
27



NURSING

EMERGENCY IV ACCESS

EASY TO USE ALERTS

WRISTBANDS

For patients with IV access placed under emergency conditions or outside the hospital.

Use at point of discharge from the emergency department - just pull red tag to activate.

Red light flashes after 18 hours prompting review of the IV access device according to hospital policy.



TTWB24H



25x250x4mm



100



IV ACCESS TUBING CHANGES - SINGLE USE LABEL ALERTS

MAKING IT EASIER TO COMPLY

LABELS

Use where you'd like to ensure compliance with IV tubing set changes.

Available in 24h, 96h and 7d timer labels.


Pull the activation tag and place on the tubing.



TTLA24H




90x25X4mm

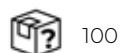

100



TTLA96H

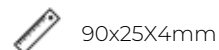


90x25X4mm



100



TTLA7D

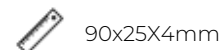

90x25X4mm



100

DRESSING CHANGES


A GENTLE REMINDER THAT IT'S TIME FOR A DRESSING CHANGE

TIMERTAG REUSABLE - 7 DAY DRESSING CHANGE ALERTS

For use in hospital wards where patients require dressing changes at 7 day intervals.

Place at the bedside or at the foot of the bed - magnetic or adhesive. Press and hold to start. Press and hold to stop.

Write on Timertag with dry-erase marker pens. Service life 12 months from date of first activation.

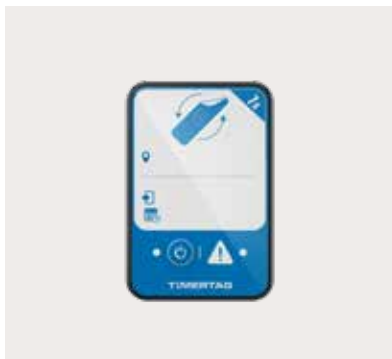


TTDC7D



63x92x8mm



20

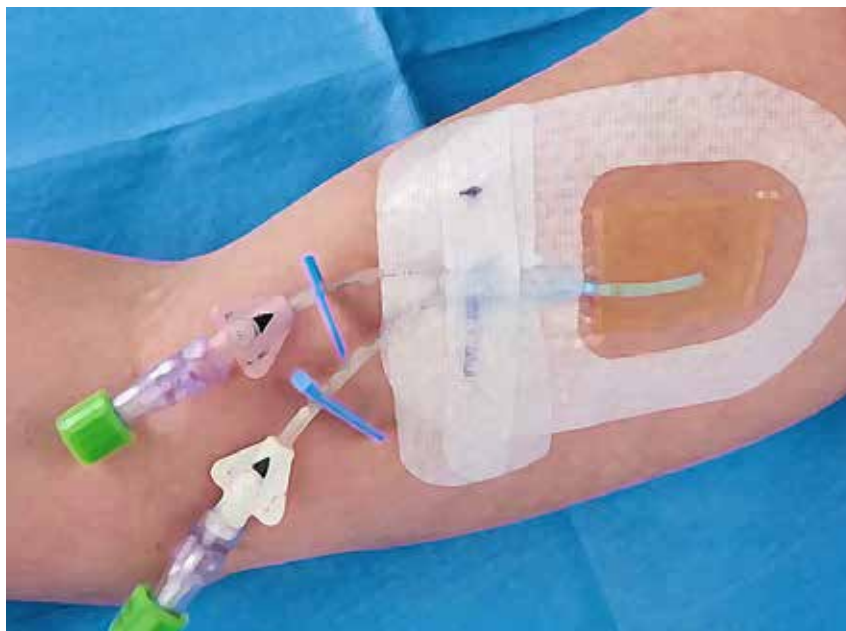

IV ACCESS REVIEWS

FOR HOSPITALS WITHOUT AN ELECTRONIC MEDICAL RECORD

TIMERTAG REUSABLE - 24H AND 8H ALERTS

For use in hospital wards where patients have IV access that requires periodic review.

Place at the bedside - magnetic or adhesive. Press and hold to start. Press and hold to stop.

Write on Timertag with dry-erase marker pens. Service life of 12 months for 24h Timertag, 6 months for 8h Timertag.



TTIV24H



63x92x8mm



20




TTIV8H



63x92x8mm



20

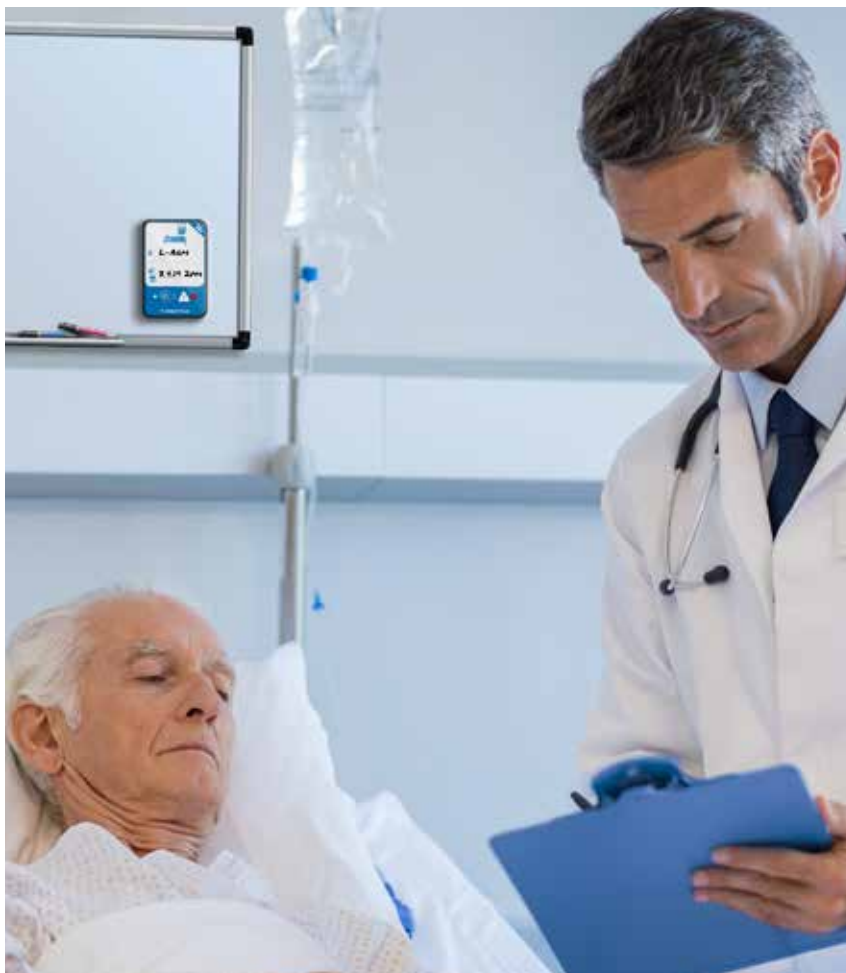

FOR HOSPITALS WITHOUT AN ELECTRONIC MEDICAL RECORD

For use in hospitals that use paper charts.

Use where patients have IV access that requires review at 72 hours.

Pull the activation tag, place in the patient chart, and tear off after use.



TTLAPIV72H



55x45X4mm



100



For use in hospital wards where patients have peripheral IV access that requires review at 72 hours.

Place at the bedside or at the foot of the bed - magnetic or adhesive. Press and hold to start. Press and hold to stop.

Write on Timertag with dry-erase marker pens. Service life of 12 months from date of first activation.

Write on Timertag with dry-erase marker pens. Service life of 12 months from date of first activation.



TTIV72H



63x92x8mm



20



ORAL CARE FOR HIGH DEPENDENCY PATIENTS

HELPING DRIVE COMPLIANCE

TIMERTAG REUSABLE - 4 HOUR ALERTS

For patients in high dependency wards that require regular oral care.

Place at the bedside or with oral care kits - magnetic or adhesive. Press and hold to start. Press and hold to stop.

Service life of 6 months from date of first activation.



TTOC4H



63x92x8mm



20



URINARY CATHETER REVIEWS - 72 HOUR

DRIVING DECISION MAKING

LABELS

Place on the urinary bag to provide an alert at point of care. Use in addition to standard daily reviews.

At 72 hours someone needs to make a decision as to whether to continue with short term catheter or not.

Pull the activation tag, place on drainage bag, and tear off after use.



TTLAUC72H



90x25x4mm



100



TIMERTAG REUSABLE - DAILY URINARY CATHETER REVIEW

For use in hospital wards without an electronic medical record, helping drive daily reviews of need for a urinary catheter.

Place at the bedside or at the foot of the bed - magnetic or adhesive. Press and hold to start. Press and hold to stop.

Write on Timertag with dry-erase marker pens. Service life of 12 months from date of first activation.



TTUC24H



63x92x8mm



20



PATIENT REPOSITIONING


A GENTLE REMINDER THAT IT'S TIME TO MOVE

TIMERTAG REUSABLE - 2 HOUR POSITION CHANGE ALERTS

For use in hospital wards where patients need repositioning every 2 hours.

Place at the bedside or at the foot of the bed - magnetic or adhesive. Press and hold to start. Press and hold to stop.

Service life of 6 months from date of first activation.

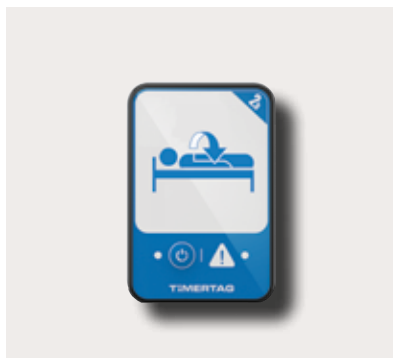

TTPC2H



63x92x8mm



20



TIMERTAG



CREATE YOUR OWN SOLUTION

MAKE NURSING EASIER YOUR WAY

TIMERTAG BLANK REUSABLE ALERTS

Use where you'd like to ensure compliance with a particular task or need.

Available in 15 different time periods from 2 hours up to 365 days - just ask.

Service life varies depending on the alert time period - anywhere from 6 months up to 5 years.



- | | |
|---|---|
|  TTBL02H |  TTBL14D |
|  TTBL04H |  TTBL28D |
|  TTBL08H |  TTBL30D |
|  TTBL12H |  TTBL90D |
|  TTBL48H |  TTBL180D |
|  TTBL24H |  TTBL365D |
|  TTBL72H |  63x92x8mm |
|  TTBL96H |  20 |
|  TTBL07D | |





 **BIONMED**

MOBILE EQUIPMENT SERVICE ALERTS

SAVE TIME AND REDUCE AUDIT RISK

BIOMEDICAL ALERTS


Use on mobile equipment that's difficult to track.

Available in 90 day, 180 day and 365 day time periods.


Pull the activation tag and stick on the equipment. "Call to action" stickers also available, or use your own.



TTBM90D

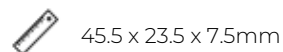


45.5 x 23.5 x 7.5mm

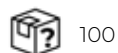


100



TTBM180D

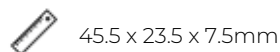


45.5 x 23.5 x 7.5mm

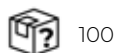

100



TTBM365D

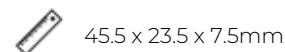

45.5 x 23.5 x 7.5mm



100



"CALL TO ACTION" STICKER





PHARMA



COMPOUNDED PHARMACEUTICAL EXPIRY ALERTS

MAKE COMPLIANCE EASIER

LABELS

For compounded pharmaceuticals sent out to wards with a 4 hour expiry time.

Pull the activation tag and place on the IV bag or delivery container.

Flashes green for four hours, then red thereafter.



TTLA04H



90x25X4mm



100



TIMERTAG





STERILE PACKS

ENSURING PACKS ARE ALWAYS READY FOR USE

LABELS

For sterile packs that are infrequently used.
If it needs a dust cover then it probably needs a Timertag.
Pull the activation tag and place on the pack.



-  TTLA365D
-  90x25X4mm
-  100



TIMERTAG REUSABLE - 365 DAY ALERTS

For sterile packs that are infrequently used.
Place on the shelf at the location where sterile pack is stored.
Press and hold to start. Press and hold to stop. Service life of 5 years from date of first activation.



-  TTBL365D
-  63x92x8mm
-  20



ENDOSCOPE STORAGE

ALERTING THAT STORAGE TIME HAS EXPIRED

LABELS

For endoscope storage.

Pull tag to activate and then loop around the endoscope. Red light flashes after set time period (7d or 21d).

Can be supplied as an adhesive label on request.



TTEN07D



250x25X4mm



100



TTEN21D



250x25X4mm



100



ENDOSCOPE STORAGE

RETROFIT YOUR DRYING CABINET WITH ELECTRONIC ALERTS

TIMERTAG REUSABLE -7 DAY ALERTS

Magnetic or adhesive attachment to drying cabinet. Press and hold to start. Press and hold to stop.

Green flashing light every 3 seconds for first 6 days, every second on final day. Red flashing light thereafter.

Service life of 3 months from first activation. Write on Timertag using dry-erase marker pen.



TTENRU07D



63x92x8mm



20



MEDICAL WASTE CONTAINERS - JCO COMPLIANCE MADE EASY

FOR CALIFORNIA, FLORIDA, PENNSYLVANIA & NEW YORK

TIMERTAG REUSABLE - 30 DAY ALERTS

Place adjacent to each medical waste container. Magnetic - metal mounting plates available.

Press and hold to start. Press and hold to stop. Red flashing light after 30 days indicates container requires changing.

Service life of 2 years from first activation. Write on Timertag using dry-erase marker pen.



TTMW30D



63x92x8mm



20



TIMERTAG

CSSD

ABOUT

INDEPENDENT RESEARCH



URINARY CATHETERS

Mitchell B et al. Reducing urinary catheter use using an electronic reminder system in hospitalized patients: A randomized stepped-wedge trial. *Infection Control & Hospital Epidemiology*, 2019; 40: 427-31.



PERIPHERAL IV ACCESS

Rutledge S et al. Clinical trial of CATH TAG: An innovative new device to reduce infections associated with indwelling medical devices. *Inf Dis Health* 2017; 22:92 - 93.



NEW INDEPENDENT RESEARCH UNDERWAY

- Urinary catheter reviews - Multicentre trial in UK, Sri Lanka, Brazil.
- IV Access compliance - Multicentre trial in Australia.
- IV Access compliance - UK Hospital.
- Oral Care - UK Hospital.
- Wristbands - US Hospital, Irish Hospital.

ABOUT

AWARDS



2016 GOOD DESIGN AWARD

"CATH TAG" PRECURSOR VERSION OF TIMERTAG LABELS



2019 DESIGN X AWARD GOLD & 2019 NEW YORK DESIGN AWARD SILVER

TIMERTAG REUSABLE ALERT



OUR

ENVIRONMENTAL STATEMENT

DISPOSAL

All Timertag products should be disposed of into battery recycling containers. The entire Timertag product is considered a “battery” for recycling purposes. We contribute to national battery recycling schemes once we reach each nation’s import thresholds.

ZERO LANDFILL

All components of Timertags (including plastic) are 100% recycled through national battery recycling schemes. Everything is repurposed and nothing goes to landfill.





ABOUT THE COMPANY

In 2013 an Infectious Diseases Physician in Australia recognised the need for an alert to remind nurses to review peripheral intravenous catheters to prevent bloodstream infections.

After discussing this with a colleague, and after failing to find an off-the-shelf solution, they decided to develop one themselves.

After successfully developing a compliance solution for peripheral intravenous catheters, the team have developed a range of new solutions for hospital care and the range of Timertag applications continues to grow.

At Timertag, we are committed to continuous development and innovation, environmental sustainability and delivering meaningful outcomes for patient care.

MEET WITH


OUR TEAM



DR GARY VERDICKT

Director

Email: gary@timertag.com

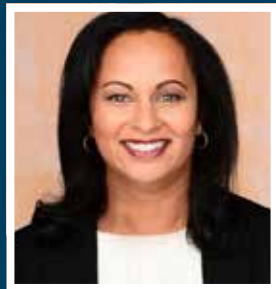

DR KUSHLANI VERDICKT

Director - Clinical Educator

Email: kushi@timertag.com



STEPHEN SONG

Business Development

Email: stephen@timertag.com



NAOMI IREDALE

Marketing Director

Email: naomi@timertag.com

PRODUCT TECHNICAL SPECIFICATIONS

LABELS, WRISTBANDS, MINI TIMERTAG

START

Pull out activation tag. Green light will flash once every second for 30 seconds. See IFU for individual products for flashing light sequence thereafter.

SHELF LIFE

5 years from date of manufacture.

PERFORMANCE

Median error (ME)=0.5%. Maximum observed error (MOE)=1.5%. TTMB90D ME=0.45 days MOE=1.35 days. TTBM180D ME=0.9 days MOE=2.7 days. TTBM365D ME=1.8 days MOE=5.4 days.

CAUTIONS

Water resistant only. Not waterproof. Do not immerse.

DISPOSAL

Dispose into battery recycling container. RoHS compliant in accordance with EU Directive 2011/65/EU.

WARNINGS

Avoid X-Rays or Magnetic Resonance Imaging (MRI). Contains a lithium manganese battery. Seek medical attention if swallowed. Keep out of reach of children. Do not burn, heat above 100°C (212°F), puncture or crush - explosion risk. Timertags are workflow aides and do not prevent disease - no health benefit claims are made.

PRODUCT TECHNICAL SPECIFICATIONS

REUSABLE TIMERTAGS

ACTIVATION

Before first use, press and hold the button until the lights flash. Product is now activated.

START

Press and hold button until green light flashes three times.

STOP

Press and hold button until red light flashes three times.

TEST

Press button for less than 1 second:

- Green blink (<1s) = Ready.
- Green light (>1s) = Active.
- Red light = Expired.
- No lights = Expired.

ALERT

Timertags 2-72h: Red light flashes every second for 2h, then once every 10s for 22h.

Timertags 7d-365d: Red light flashes every second for 8h, then once every 10s for 14d.

SHELF LIFE

Product must be activated within 2 years from date of manufacture.

SERVICE LIFE

Timertags work an unlimited number of times from first activation for the following time periods:

- <12h Timertags - 6 months.
- 24h-7d Timertags - 12 months.
- 30d Timertags - 2 years.
- 90d Timertags - 3 years.
- 1y Timertags - 5 years.

If Timertag is not activated within 2 years of manufacturing date Timertag will not function.

CAUTIONS

Water resistant only. Not waterproof. Do not immerse. Only use dry-erase (whiteboard) pens on Timertags.

DISPOSAL

Dispose into battery recycling container. RoHS compliant in accordance with EU Directive 2011/65/EU.

WARNINGS

Avoid X-Rays or Magnetic Resonance Imaging (MRI). Contains a lithium manganese battery. Seek medical attention if swallowed. Keep out of reach of children. Do not burn, heat above 100°C (212°F), puncture or crush - explosion risk. Timertags are workflow aides and do not prevent disease - no health benefit claims are made.



TIMERTAG

PO BOX 1187, SUTHERLAND NSW 1499 AUSTRALIA
+61 2 9521 8880 info@timertag.com timertag.com